

Ocean SR Series: SR6 Spectrometer



Great Versatility, High Sensitivity

The **SR6** is a compact, UV-Vis-NIR spectrometer anchored by a 2048-element CCD-array detector and powerful electronics to provide high spectral response, high optical resolution and excellent thermal wavelength stability in a single package. SR6 is an ideal choice ranging from

lab settings with controlled environments to industrial locations where temperature fluctuations can affect instrument performance.

US +1 727-733-2447


EUROPE +31 26-3190500

ASIA +86 21-6295-6600

info@oceanoptics.com • www.oceanoptics.com



OceanOptics

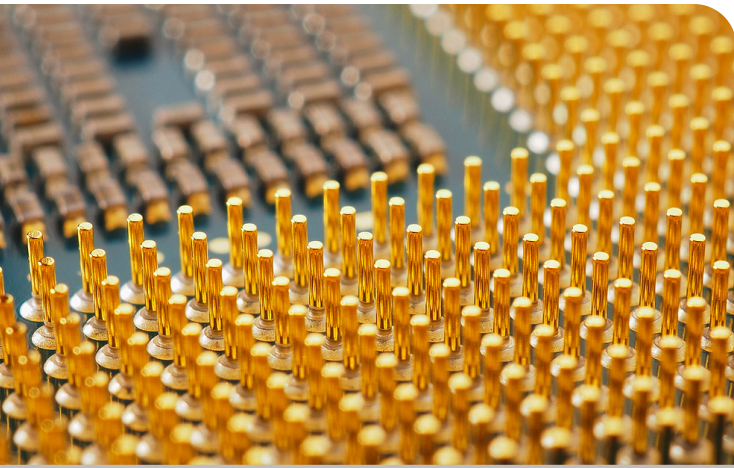

At a Glance

Wavelength Range: ~180 - 1100 nm

Optical Resolution (w/25 μ m slit): 0.60-2.03 nm

Dynamic Range: 12000:1

Integration Time: 7.2 ms - 5 s

Signal-to-Noise Ratio (max. per second w/ High Speed Averaging Mode): 3500:1

Signal-to-Noise Ratio (sigle scan): 400:1

Scan Rate: 74 scans / second

Thermal Stability: 0.02 nm / ° C

Interfaces: USB, RS-232, Ethernet

Connector: USB Type-C, Samtec TFM-108-02-L-DH, RJ45

Temperature (storage): -30 ° C-70 ° C

Temperature (operation): 0° C-55° C

A/D Resolution: 16

High Speed Average Mode: Yes

Dimensions: 88.1 mm x 63.5 mm x 31.4 mm


Weight: 254 g

Versatility of the SR6

SR6 is a multipurpose spectrometer, valued as both a one-off setup for the lab or as a customized system for high-volume industrial and OEM applications. The spectrometer comes in models covering wavelength ranges within ~180-1100 nm, offers entrance slit options in widths from 5 μ m to 200 μ m, and couples to light sources, optical fibers and sampling optics to optimize configurations for various applications.

SR6 spectrometers operate with OceanView spectroscopy software and include OceanDirect, a powerful, cross-platform Software Developers Kit with an Application Programming Interface. With its library of functions, OceanDirect allows users to write custom software solutions for their spectrometer, optimize spectrometer performance and access critical data for analysis. OceanDirect also enables High Speed Averaging Mode, a hardware-accelerated signal averaging tool that markedly enhances spectrometer signal to noise ratio (SNR) per unit time. This promotes higher quality spectra and more accurate results.

Additional Functionality with Ocean SR6

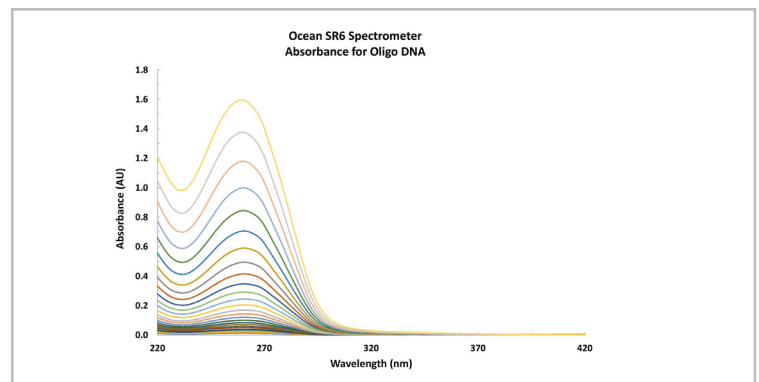

The SR6 comes in multiple UV models. Here's an example showing absorbance of oligo DNA concentrations.

SR6 has trigger mode options that enable actions such as synchronizing spectral acquisition to an external event (e.g., pulsing of a lamp) or timing spectral acquisition to meet certain sampling conditions. Triggering provides accurate timing and synchronization between the SR6 spectrometer and other devices, adding another layer of versatility to your spectral measurements.



For more information on the SR6, please contact an Ocean Optics Application Scientist today.